



UniversalLine2

Technology for people

Functionality and versatility.

Relaxed use with optimal comfort.





UniversalLine2

Technology for people

The UniversalLine2 is suitable for outpatient treatments, dialysis, oncology, blood and plasma donation.

All aspects of clinical practice are taken into account.

All segments can be adjusted for optimal relaxation.

The UniversalLine2 also has an easily adjustable, soft head cushion. The manually adjustable footrest further ensures a comfortable sitting position.


QUALITY
MADE IN
GERMANY



The UniversalLine2 is equipped with 4 motors that enable quick and smooth adjustment of the back, seat and leg sections as well as height adjustment.



The illuminated handset impresses with its instinctive operation.

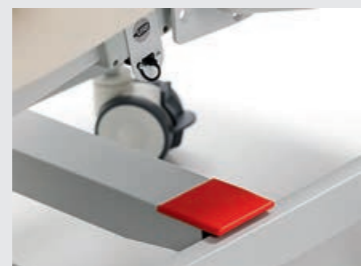

The large foot switch allows the head-down position to be achieved rapidly and with ease.



All motors are designed for minimal maintenance and ease of servicing.

The comfortable pillow can be easily adjusted in height.



Give your surroundings more colour

The high-quality design of the therapy table with your choice of upholstery color makes it an ideal addition to any hospital, clinic or doctor's surgery.

A wide range of faux leather colors also helps you to create a pleasant and relaxing atmosphere in your environment.

The pleasant feel and excellent material properties ensure long-lasting, easy-care and comfortable use.

The adjustable footrest helps to minimize discomfort.

Equipotential bonding increases patient safety.

Hygiene in perfection

Smooth surfaces and a three-part lying surface ensure quick and easy cleaning - even in the spaces in between. The absence of Velcro as a fastening method also significantly enhances hygiene conditions.

The upholstery is also permanently protected against staining with a PVC protective cover. The protective cover is available for the leg section, the armrest and the headrest cushion.



Innovative overload protection

Thanks to the innovative overload protection system, the armrest springs back when the load exceeds 25 kg, thus increasing safety for patients when getting into the chair. At the same time, the overload protection prevents damage reducing repair costs and downtime. It can be moved to the desired position with just one lever position and rotate completely so that staff can use it instinctively.



Technical data

Dimensions

Width incl. armrests	90 cm
Width of the seat	58 cm
Overall length	203 – 230 cm
Armrests	52 x 15,5 cm
Total weight	94 kg

Adjustment options

Backrest angle	75° – 0°
Seat section angle	30° – 5°
Leg section angle	5° – 30°
Footrest	28 cm
Height adjustment	54 – 80 cm

Safety

Safe working load	240 kg
-------------------	--------



Manufacturer:
MEDICAL GmbH Medizintechnik
Ravensberger Straße 65
32278 Kirchlingern, Germany



Distributed by:
MEDICAL Sales GmbH & Co. KG
Wilhelmstraße 14 A
61381 Friedrichsdorf, Germany

Individual accessories



Reading lamp



Handle on backrest



IV-pole



Cover leg segment



Console table



Paper roll holder



USB port



Rechargeable battery



2nd hand switch

Note: The product shown in this brochure contains optional accessories that are not included in the standard version.